

New designs for more dynamic outcome Carbon pre-preg ankle-foot-orthosis


Streifeneder
ortho.production GmbH
Moosfeldstrasse 10
82275 Emmering
Germany
T +49 8141 6106-0
F +49 8141 6106-70
export@streifeneder.de
www.streifeneder.com/op



300W125/E 2018-04



The PeroSupport.tec

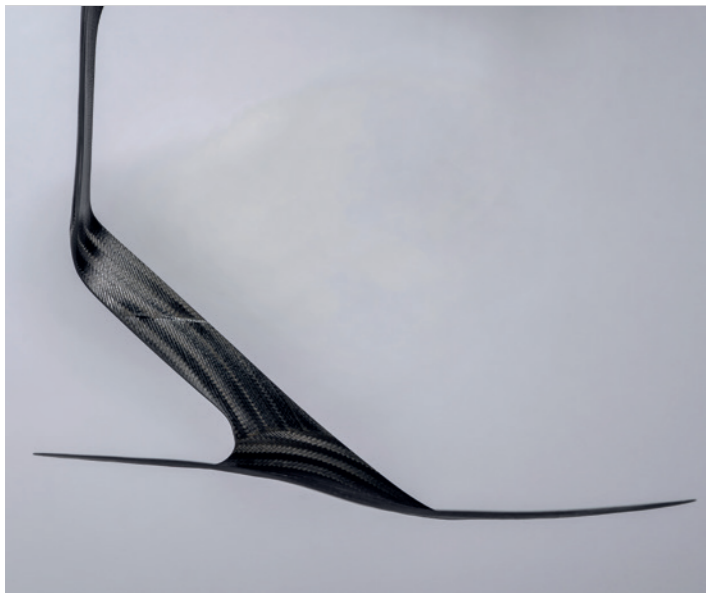

Flexible lightweight

The new PeroSupport.tec is made of high-quality carbon pre-preg what is extremely stable and lightweight. The material stores the energy gained in the stance phase and gives it back to the user in the swing phase. The forefoot and heel areas are flexible to allow the patient the most physiological and dynamic gait possible. The calf section and metatarsal area are load-bearing elements that actively support foot lifting.



Perfect base

The reduced overall height of the PeroSupport.tec and the improved sole ensure maximum wearing comfort: The lateral design (heel-pitch) of the new model is much more pronounced in order to improve the adaptation to the patient's shoe. The sole shape itself has become narrower and can now be cut to the marked outer edge with sturdy scissors.



Donning with a single click

Simple donning and doffing is particularly important for drop foot orthoses. Many users are limited in their fine motor abilities due to hemiparesis (apoplex/stroke). The new, high-quality magnetic closure of the Fidlock brand perfectly meets these requirements. The closure engages by itself with a soft click and the patient immediately knows whether the orthosis is closed securely. Opening is also quick and uncomplicated thanks to the magnet.

New:

The Fidlock® Magnetic Closure 44P20/25 is also available seperately as flexible spare part.

Indications

- peroneal paresis
- weakness of dorsiflexion of the foot
- inflammation of the ankle joint
- postoperative treatment after achilles tendon rupture

Good to know:

The sole of the PeroSupport.pro ankle-foot-orthosis can now also be cut with sturdy scissors. It offers additional rotation and calf stability.

PeroSupport.tec

Size	Side	Shoe size	Height	Item-No.
XS	Left	34 – 36	32.5 cm	197R4/L0
S	Left	36 – 39	33.5 cm	197R4/L1
M	Left	39 – 42	34.5 cm	197R4/L2
L	Left	42 – 45	35.5 cm	197R4/L3
XL	Left	45 – 48	37.5 cm	197R4/L4
XS	Right	34 – 36	32.5 cm	197R4/R0
S	Right	36 – 39	33.5 cm	197R4/R1
M	Right	39 – 42	34.5 cm	197R4/R2
L	Right	42 – 45	35.5 cm	197R4/R3
XL	Right	45 – 48	37.5 cm	197R4/R4