

One step ahead.
The AeroWalker and AeroWalker.short



The rigid shell construction of both pneumatic Ankle-Foot-Orthoses AeroWalker 68R1 (long version) and AeroWalker.short 69R1 (short version) provides controlled immobilization in 90° functional position of the lower leg, ankle and foot.



Lightweight synthetic shells and functional design

With their extremely low weight, the AeroWalker as well as the AeroWalker.short make daily living of your patient much easier and therefore increase their life quality during therapy. Because of the slotted construction of the shells, the orthoses are highly breathable and comfortable. Due to the improved rolling motion, walking with the AeroWalker and AeroWalker.short is made easier



Breathable insole

Every AeroWalker and AeroWalker.short is now equipped with a new, breathable insole. Small perforations allow the air to optimally circulate and prevent possible heat accumulation. In case humidity should still occur, it will escape directly instead of accumulating underneath the foot sole. The user's injured foot can therefore recover under optimal and comfortable temperature conditions.

More comfort, better fit: the anatomically shaped insole blank can be individually adjusted - for example by grinding - and effectively supports the transverse and longitudinal arch of the foot. It is available in different sizes and comes with a black microfleece cover fabric.



Size	Length	Width	Item-No.
35 – 36	29.5 cm	11.5 cm	62R15/35...
37 – 38	30.5 cm	12.0 cm	62R15/37...
39 – 40	31.5 cm	12.5 cm	62R15/39...
41 – 42	33.0 cm	13.0 cm	62R15/41...
43 – 44	34.0 cm	13.5 cm	62R15/43...
45 – 46	35.5 cm	14.0 cm	62R15/45...
47 – 48	36.5 cm	14.5 cm	62R15/47...

Ordering example: "62R15/35" + "L" (left side) or "R" (right side)



Simple and comfortable use


Integrated in the inner shoe, an inflatable air chamber completely embraces the foot from ankle joint to midfoot. This provides a perfectly stabilized position of foot and lower leg. The simple and comfortable use of the air pump allows individual inflation of the air chambers, and therefore provides optimal wearing comfort. A rounded, slightly bent to the outside, slotted edge in the calf area helps to prevent pressure marks.



Stable, safe and comfortable for daily use

The stable and rigid sole of the AeroWalker and the AeroWalker.short provides immobilization of foot- and toe joints and with its rocker-bottom sole shape, supports the natural rolling motion. The adjustable and removable toe protection serves to avoid concussion damages in the front area. The profiled anti-slip sole grants additional safety during walking. The insole is optimally cushioned for increased shock absorption during heel strike.

Position number: 23.06.01.0063 (AeroWalker),
23.06.01.0064 (AeroWalker.short)

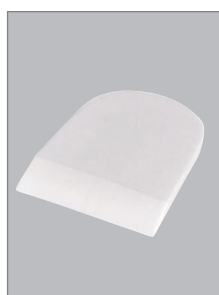

With additional application of the heel wedges (62R13 / 5-pc. set), the AeroWalker and AeroWalker.short can also be used for load relief after injuries of the achilles tendon.

Size	Shoe size	AeroWalker Item-no.	AeroWalker.short Item-no.
S	35 – 37.5	68R1/S	69R1/S
M	38 – 42.5	68R1/M	69R1/M
L	43 – 45.5	68R1/L	69R1/L
X	> 46	68R1/XL	69R1/XL

Inner shoe 68R11 (for AeroWalker) or 69R11 (for AeroWalker.short) and two padded socks are included.

Streifeneder
ortho.production GmbH
Moosfeldstraße 10
82275 Emmering
Germany
T +49 8141 6106-0
F +49 8141 6106-50
office@streifeneder.de
www.streifeneder.de/op



300W127/E 2021-07